



Main Unit

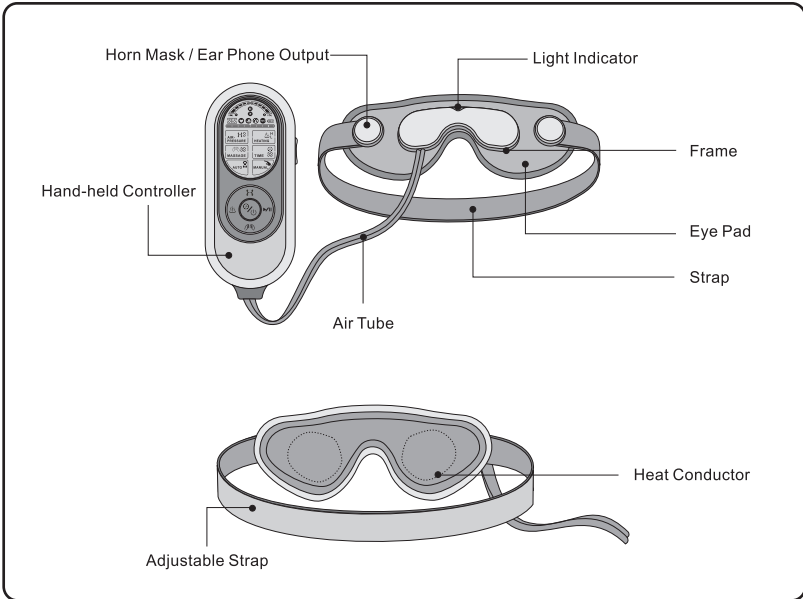

Accessories

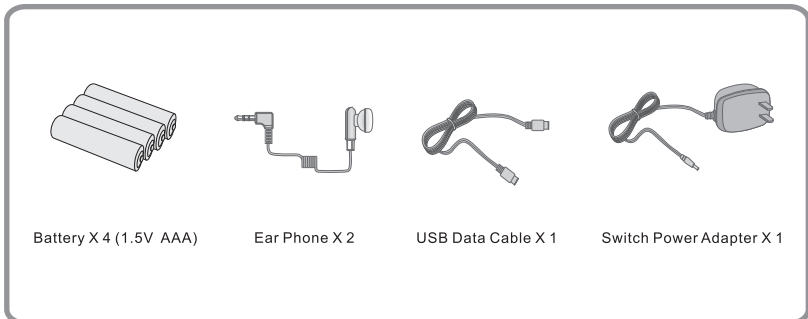


Power Supply


- 4X1.5V AAA battery (DC 6V) or AC power adapter (Output DC 6V)
- **Before using the data cable, please make sure that the dry cell batteries or the power supply is disconnected.**

Inserting the batteries


❶ Press the back cover and slide it upwards, according to the arrow

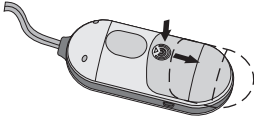


❷ Insert the batteries correctly according to the polarities

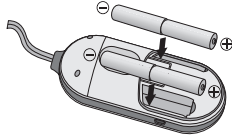

❸ Put the battery cover back

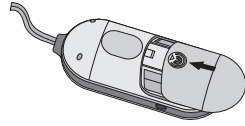


Connecting it to the power supply

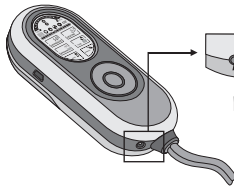


❶ Connect the plug from the AC power adapter to the base of the controller.

❷ Connect the AC power adapter to an AC socket

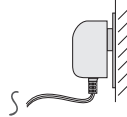

Function Keys

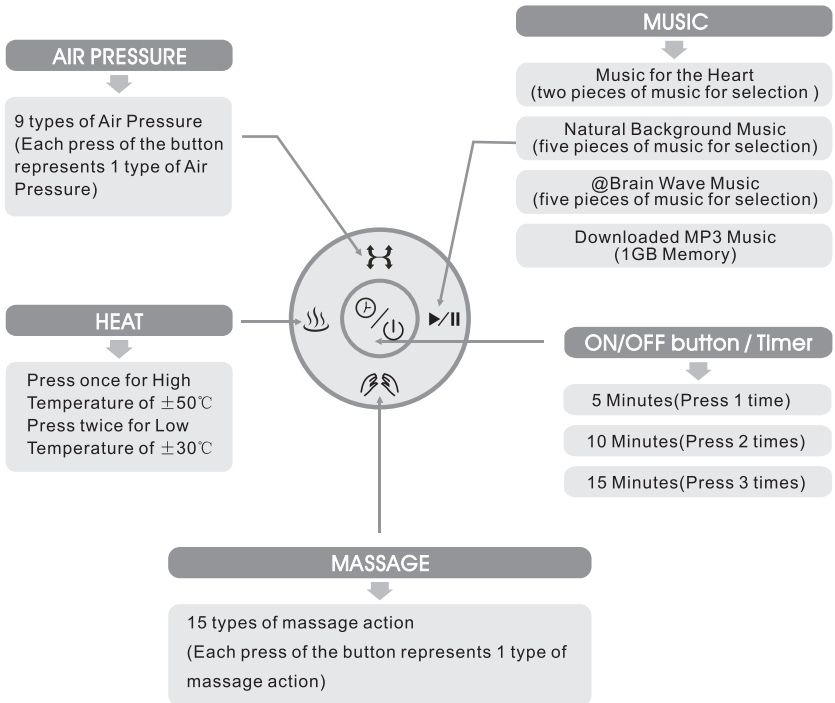

LCD Screen

Symbols	Function	Description
	Air Pressure	Types 1-9
	Heat	High Temperature of $\pm 50^{\circ}\text{C}$ Low Temperature of $\pm 30^{\circ}\text{C}$
	Massage	Types 1-15
	Timer	5, 10 or 15 minutes
	Power	Low power or weak batteries when it flickers
	Play/Pause	To start/interrupt playback of music
	Volume	Volume Adjustment Display
	Music N	Number of songs being displayed
	Music for the Heart	2 pieces of music for selection
	Natural Background Music	5 pieces of music for selection
	@Brainwave Music	5 pieces of music for selection
	MP3	Downloaded MP3 Music
	Music Progress Bar	Progression Display

● Note:


1. Only the icon of the type of song being played is shown, while the other 3 will not be displayed.
2. When either icon is selected, the sound track chosen will continue to repeat itself in a loop. To skip to the next music sound track, please press buttons to switch to the previous or to the next song, However, when icon is chosen, the music sound track will automatically skip to the next song upon finished playing the current song.

Operating Instructions

Automatic Mode

- Press the **ON/OFF/TIMER** button and the machine will enter the 3 different programmes.
- It operates for 5 minutes under each combination mode, after which the order of combinations will alternate automatically .

Programmes	Time
AIR PRESSURE / MASSAGE / MUSIC	5 minutes
MASSAGE / HEAT / MUSIC	5 minutes
HEAT / AIR PRESSURE / MUSIC	5 minutes

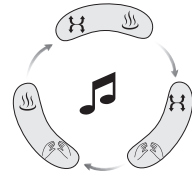

Manual Mode

Once you press the ON/OFF button to start the machine, the machine will automatically enter the Automatic Mode. To disable Automatic Mode and enter the Manual Mode, simply press and hold that particular Therapy button for around 2 seconds to disable that function.

Example:

The massager automatically enters the first Automatic programme: Air Pressure / Massage / Music.

If the user only wants Air Pressure Therapy and not Massage Therapy or Music Therapy, simply press and hold the Massage Therapy button and Music Therapy button each for around 2 seconds to disable that function. And the machine will only operate Air Pressure Therapy.


Similarly, if the user wants the Massage function but not the Air Pressure Therapy or Music Therapy, simply press and hold the Air Pressure button and the Music Therapy button each for around 2 seconds to disable that function. And the machine will only operate Massage Therapy. Perform likewise if the user wants to select the Music Therapy only.

Operating Instructions


Downloading MP3 Music

Connecting it to the data cable (For MP3 music purpose)


- 1 Please insert USB data cable into the USB socket of the controller

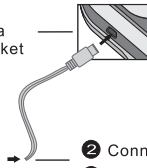


- 2 Connect to the data source
- 3 You can start downloading your desired MP3 songs

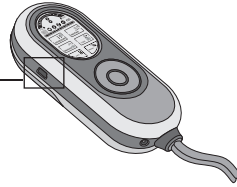

MP3 Function Keys

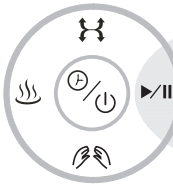

1. Press and hold for around 2 seconds to enable or disable Mp3 function;
2. Press the button to switch between Play & Pause.

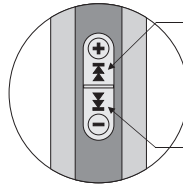


- Press this button to skip to the next song
- Press & hold for around 2 seconds to increase the volume.
- Press this button to play the previous song
- Press & hold for around 2 seconds to reduce the volume.

● **Remarks:**

You can still listen and enjoy the downloaded MP3 music even when you do not want any heat / air pressure / massage actions from the machine.

Explanation of MP3 icons

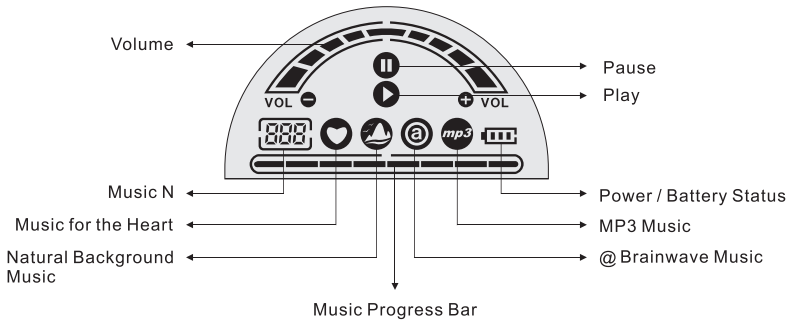

Troubleshooting Guide

Problem	Possible causes	Solution	Remarks
The eye massager fails to operate	<ol style="list-style-type: none">① Weak batteries② Batteries installed incorrectly③ The batteries are not firmly affixed in the battery compartment	<ol style="list-style-type: none">① Replace with new batteries & check the cable connection② Ensure the polarities of the batteries are correct③ Ensure the batteries are firmly affixed in the battery compartment	Contact your local authorized dealer if problem persists
Stimulation is too weak	Weak batteries	Replace with new batteries	
LCD displays fuzzy or no image	Weak batteries	Replace with new batteries	

Specifications

Product Size

- Eye massager : 230x92x105mm
- Hand -held controller : 138x60x42 mm
- Net weight of main unit: 360 g

Operating Environment

- Temperature : 5°C~ 40°C
- Humidity : ≤80%
- Atmospheric pressure : 86kPa~106kPa
- AC Power Adapter : AC100-240V 50/60Hz / DC6V 0.70A
- Battery : 4X1.5V AAA